

Product Specification

FEATURES AND BENEFITS

- Accepts Chemetron®, Diamond, Puritan Bennett® and DISS specific adapters
- Pin indexed to prevent interchangeability of gas service
- Each outlet is cleaned for medical gas service and pressure tested
- Less than 3 psi (21 kPa) pressure drop through the outlet at 120 L/min and 50 psi (345 kPa) inlet pressure
- Gas specific back bodies can accept either Quick Connect or DISS front identification bodies
- Outlets can be adjusted up to 3/4" (19mm) in mounting plate thickness
- Adapters are available from Ohio Medical Corporation

SPECIFICATION

The medical gas ceiling column outlet shall be Ohio Medical® series. The outlet shall be UL® Listed, NFPA® 99 compliant, cleaned for medical gas service and be pressure tested. Each outlet shall have less than 3 psi (21 kPa) pressure drop through the outlet @ 120 L/min, and 50 psi (345 kPa) inlet pressure. For outlets providing positive pressure gas, the outlet shall be equipped with a primary and secondary check valve, the secondary check valve shall be rated for 200 psi (1,379 kPa) allowing the primary check valve to be removed for service without isolating the entire zone.

The ceiling column outlets shall have a gas specific back body which can accommodate up to a 3/4" (19mm) plate. Each back body shall be equipped with a 6-1/2" (165mm) length type "K" copper pipe stub which is brazed to the outlet body. The outside diameter of the copper pipe stub shall be 1/2" (12.7mm). The body shall be 1-5/16" (33mm) diameter one piece brass construction.

Outlet bodies shall be gas specific by means of a gas specific indexing pin arrangement. The indexing pins shall be arranged in such a way as to connect the latch-valve assembly only with the specific matching gas back body, preventing interchangeability of gas services.

The latch-valve assembly, which by means of color coding and wording shall identify the specific medical gas service provided by the outlet and shall accept gas specific adapters.

Equipment is manufactured in an ISO® 9001 registered facility.



Diamond Compatible Console Outlet




Chemetron Compatible Console Outlet



D.I.S.S. Compatible Console Outlet



Rough-in assemblies to be mounted to column. Piping can not support weight of outlet assembly.

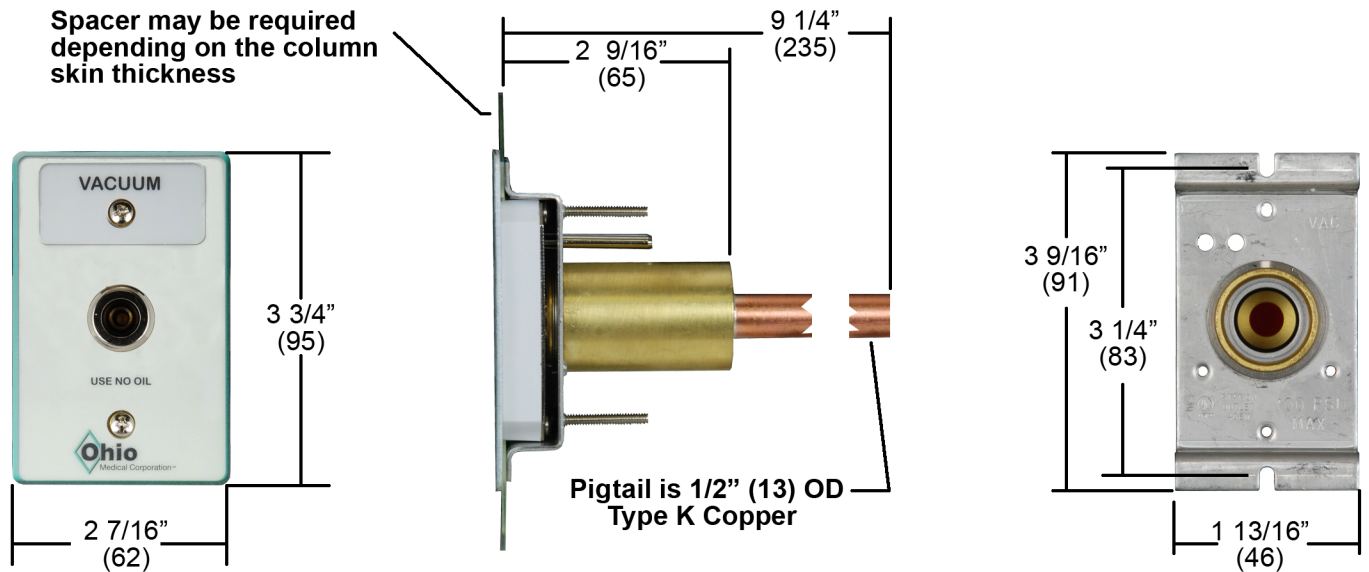

Dimensions in Inches (mm)

OUTLET DESCRIPTION	DISS CEILING COLUMN	DIAMOND COMPATIBLE CEILING COLUMN	CHEMETRON COMPATIBLE CEILING COLUMN	PURITAN BENNETT COMPATIBLE CEILING COLUMN
Oxygen	261000-4	261010-4	261020-4	261040-2
Vacuum	261000-8	261010-8	261020-8	261040-6
Air	261000-16	261010-16	261020-16	261040-14
Nitrous Oxide	261000-12	261010-12	261020-12	261040-10
WAGD	261000-24	261010-20	261020-20	261040-18
Nitrogen	261000-20	N/A	N/A	N/A
Carbon Dioxide	261000-28	N/A	N/A	N/A
Instrument Air	261000-32	N/A	N/A	N/A

© 2015 Ohio Medical Corporation.
 This document contains information that is proprietary and confidential to Ohio Medical Corporation. Use of this information is under license from Ohio Medical Corporation. Any use other than that authorized by Ohio Medical Corporation is prohibited.
 Ohio Medical Corporation, Ohio Medical and the Ohio Medical Corporation logo are registered trademarks of Ohio Medical Corporation.
 Chemetron is a registered trademark of Allied Healthcare Products, Inc., Puritan Bennett is a registered trademark of Nellcor Puritan Bennett LLC.
 UL and the UL logo are registered trademarks of UL LLC Limited Liability Company. NFPA is a registered trademark of National Fire Protection Association, Inc.
 ISO is a registered trademark of International Organization for Standardization